

NEW

lighted **BIG CASE back tables**



Now Features:
LED lights and electric lift system

Big Case Back Tables

Tables are designed to support a full working weight capacity of 550 lbs (249.5kg) distributed evenly.

Top Tier
150lbs (68kg)

Main Surface
200lbs (91kg)

Bottom Storage
200lbs (91kg)

Standard
427 4' table (1.5m)
428 5' table (1.5m)
429 6' table(1.8m)

Adjustable Electric/Lighted
427-ADJ-E 4' table (1.5m)
428-ADJ-E 5' table (1.5m)
429-ADJ-E 6' table(1.8m)

Lighted
427-L 4' table (1.5m)
428-L 5' table (1.5m)
429-L 6' table(1.8m)

Top Tier Height Adjustment

15" (38cm)

12" (30cm)



Main Surface
Height Adjustment



▲ Height Adjustment ▼


The UP/DOWN buttons for height control are located to the right of the power supply. From the lowest position to the maximum height takes approximately 55 seconds.





l i g h t e d

- LED light strip is fully enclosed and water resistant
- Easy to clean
- 6300K cool white light for easy cleaning and durability.

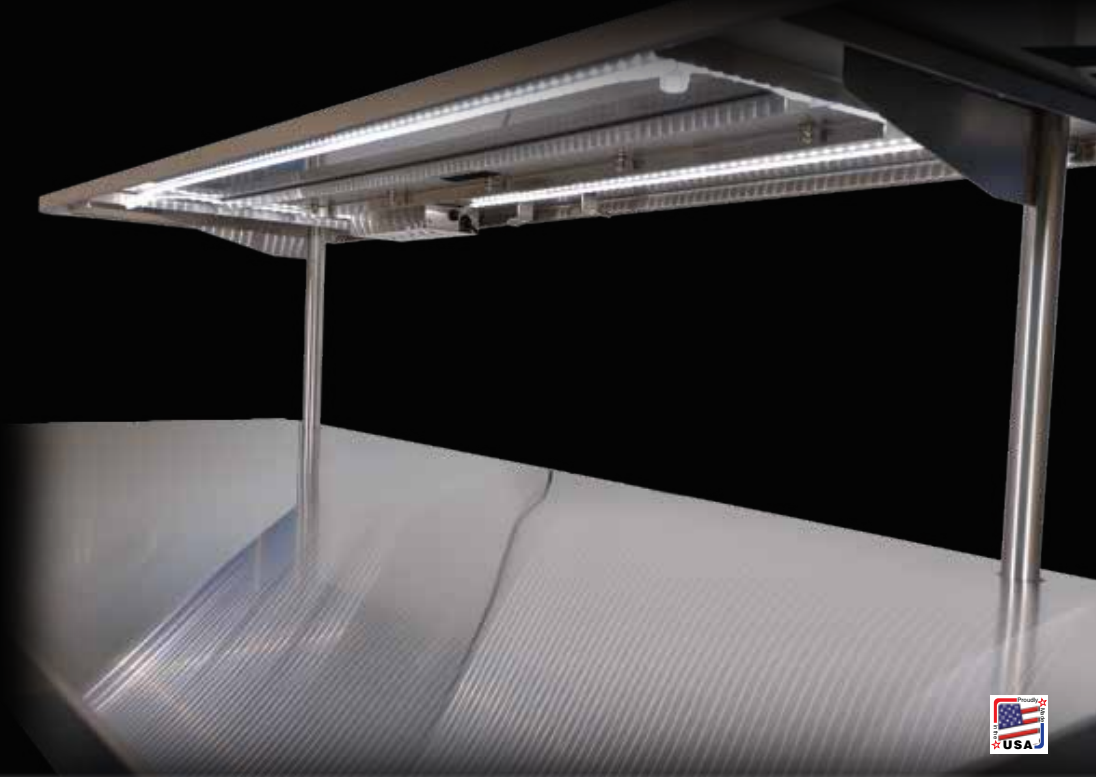

BIG CASE back tables



Standards

For patient and staff safety, the Big Case Back Table meets the highest hospital safety standards. IEC 60601-1 and UL approved E 503 186 approved.



Power Cord

Hospital grade power cord can be wound around storage loops when not in use. Cord length measures 9' 8" (3m).




Power Supply

Located at the back of the table are two independent ON/OFF toggle switches. One to control the entire table unit and a second toggle for the LED lights.



Casters

Electric adjustable back tables feature four Total Lock Casters measuring 4" (10cm) diameter.

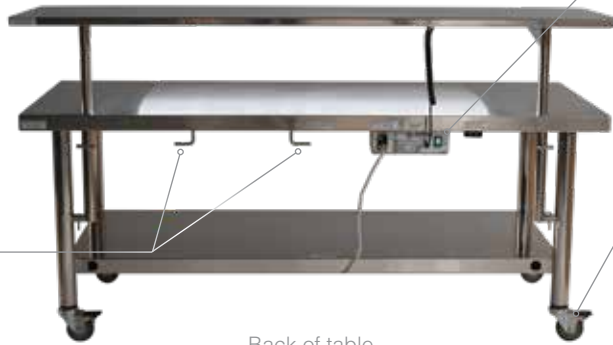

Back of table



Trusted Fit Disposable Drapes

Unique absorbent drape easily covers both the upper and main surfaces of the table. The clear plastic backdrop protects the sterile field while allowing light through to the main table surface

418-HD-S 4' table (1.5m) 10 per case

419-HD-S 5' table (1.5m) 10 per case

420-HD-S 6' table (1.8m) 10 per case



Upper tier shown in storage position

