

Two-Tier Table Drapes (Protected by U.S. Patent No. 6,019,102)

- Unique single piece disposable drapes:
Heavy Duty: Three layer system including extra absorbent center layer - 28 mil thick.
- One piece drape easily covers both upper and main table surfaces.
- The clear plastic backdrop protects the sterile field while allowing light through to the main table surface.



TABLE & DRAPE CROSS REFERENCE CHART			Perfect-Fit Drapes for OR Specific's Big Case Back Tables				
Part No.	Description	Dimensions	Heavy-Duty				Per Case
			418-HDE-S	419-HDE-S	419-HDE-NS	420-HDE-S	
427	122 cm Back Table with removable 55 cm D second tier shelf	122 L x 76 P x 84 H	•				10
427-ADJ	122 cm Back Table with removable 55 cm D second tier shelf (Height adjustment from 74 cm to 100 cm)	122 L x 76 P x 74 H	•				10
428	152 cm Back Table with removable 55 cm D second tier shelf	152 L x 76 P x 84 H		•	•		10
428-ADJ	152 cm Back Table with removable 55 cm D second tier shelf (Height adjustment from 74 cm to 100 cm)	152 L x 76 P x 74 H		•	•		10
429	183 cm Back Table with removable 55 cm D second tier shelf	183 L x 76 P x 84 H				•	10
429-ADJ	183 cm Back Table with removable 55 cm D second tier shelf (Height adjustment from 74 cm to 100 cm)	183 L x 76 P x 74 H				•	10



Draping Video DVD Demo

Request a copy of the features, benefits and draping instructions for the Big Case Back Table. Download available at www.orspecific.eu.



Efficiency Report

Request a copy of the time saving benefits for the Big Case Back Table from your local rep. PDF download available at www.orspecific.eu.



30 Day Evaluation

To learn more about our table and drape purchasing program, please contact OR Specific's customer service department at 800.937.7949.



MMORS-50-EU
01/20

BIG CASE back tables

Specifically designed for a variety cases



Maximize Space

Upper shelf creates additional space without compromising visibility



Visibility

Easily arrange & organize instrument trays without the need for stacking



Time

Improve setup efficiency & workflow during large surgical cases



Drape

Create a sterile field with our patented one-piece drape.

Standard or heavy-duty
Protected by U.S. Patent No. 6,019,102



standard BIG CASE Back Tables

The OR Specific standard Big Case Back Table is ideal in supporting any surgical procedure that is implant and/or instrument intensive, especially when space is limited.

BIG CASE back tables

adjustable BIG CASE Back Tables

OR Specific's adjustable back tables accommodate variable staff heights ranging from 149 cm to 198 cm. The main surface of the table may be adjusted intraoperatively from a minimum height of 74 cm to a maximum height of 100 cm.

STANDARD TABLE SPECIFICATIONS

The standard Big Case Back Table has a fixed table height of 84 cm (floor to main surface). When upper tier is in the storage position, the table can be used as a single tier back table.

- **Upper Shelf Set Positions**
30 cm above main surface
38 cm above main surface
- **Weight Capacity**
68 kg upper shelf
90 kg main surface
90 kg storage shelf
— weight uniformly distributed
- **18 gauge stainless steel**
- **Four 9 cm casters**
- **Two front locking casters**



427 Big Case 122 Standard

Table Dimensions (cm): 122 L x 76 D x 84 H


Shelf Dimensions (cm): 122 L x 55 D



428 Big Case 152 Standard

Table Dimensions (cm): 122 L x 76 D x 84 H

Shelf Dimensions (cm): 152 L x 55 D

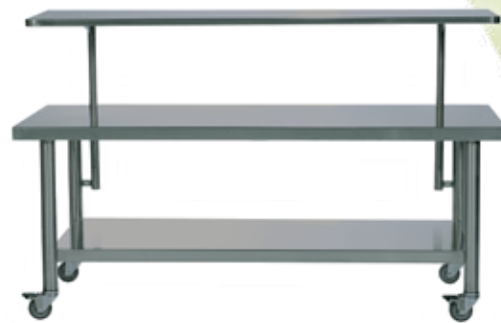

429 Big Case 183 Standard

Table Dimensions (cm): 183 L x 76 D x 84 H

Shelf Dimensions (cm): 183 L x 55 D



427-ADJ Big Case 122 Adjustable

Table Dimensions (cm): 122 L x 76 D x 74 H
Adjustable (cm): 74 LOW to 100 HIGH

Shelf Dimensions (cm): 122 L x 55 D



428-ADJ Big Case 152 Adjustable

Table Dimensions (cm): 152 L x 76 D x 74 H
Adjustable (cm): 74 LOW to 100 HIGH

Shelf Dimensions (cm): 152 L x 55 D



429-ADJ Big Case 183 Adjustable

Table Dimensions (cm): 183 L x 76 D x 74 H
Adjustable (cm): 74 LOW to 100 HIGH

Shelf Dimensions (cm): 183 L x 55 D

ADJUSTABLE TABLE SPECIFICATIONS

Adjustable Big Case Back Table models may be raised or lowered intraoperatively through the use of a hand crank.

- **Adjustable Main Surface**
74 cm minimum height
100 cm maximum height
- **Upper Shelf Set Positions**
30 cm above main surface
38 cm above main surface
- **Weight Capacity**
68 kg upper shelf
90 kg main surface
90 kg storage shelf
— weight uniformly distributed
- **18 gauge stainless steel**
- **Four 9 cm casters**
- **Two front locking casters**



Two locking casters on front legs of table hold table stationary.




Table shown with upper tier in storage position




Secure the upper tier with T-bolts once raised and positioned.




Raise or lower the table by using the hand crank at back of table.

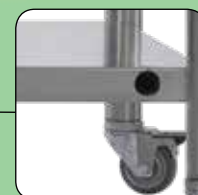

Post hole for storing upper tier when not in use.

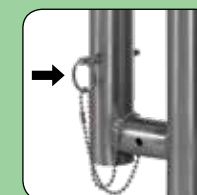

Safety locking pin used for securing upper tier shelf.

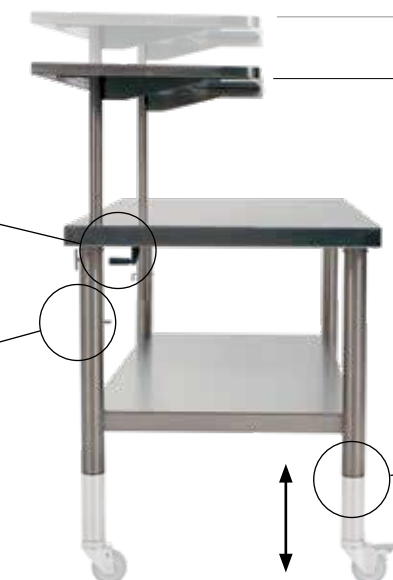

Set at 38 cm above main surface

Set at 30 cm above main surface

100 cm max. height

74 cm min. height

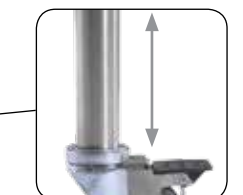


Table has a 23 cm vertical travel to reach a maximum height of 100 cm.